

# REFRIGERATOR & FREEZER

# ACCESSORIES

(See Price Sheet for Additional Accessories)

## Digital Temperature Displays and Alarms

### TDA-120 RC

120 VAC operated freestanding digital temperature display with high and low alarm (audible and visual), remote contacts, mute button and alarm delay. Can be surface mounted to most refrigerator and freezer units or nearby walls or work surfaces. Includes 10' power supply cord and 10' of sensor cable to accommodate temperature probe installation within refrigerator or freezer. LCD display for clear reading of current temperatures. Requires simple field installation.


### TDA-120 RC-4

Same as TDA-120RC but with the capability of monitoring up to 4 channels (comes standard with 2 probes)


### TDA-B

Battery operated freestanding traceable memory monitoring thermometer with high and low alarm, perfect for monitoring temperatures in refrigerators and freezers. Minimum and maximum temperature memory. Triple display simultaneously shows high and low and current temperatures. Temperature probe is held in glycol filled bottle for more realistic sample readings. Supplied with a 10' micro-cable to permit easy installation through refrigeration and freezer gaskets. Can be surface mounted directly to refrigerator or freezer. Includes 1 AA alkaline battery for continuous, 1 year monitoring. Traceable to NIST for security, certificate included. Requires simple field installation.


TDA-120 RC

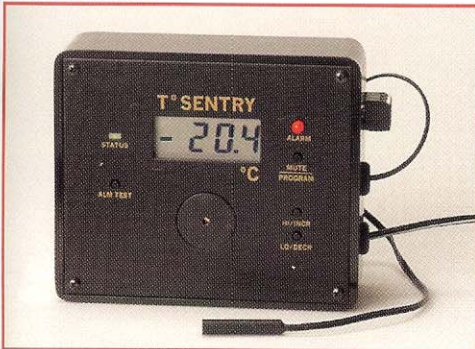

TDA-120 RC-4

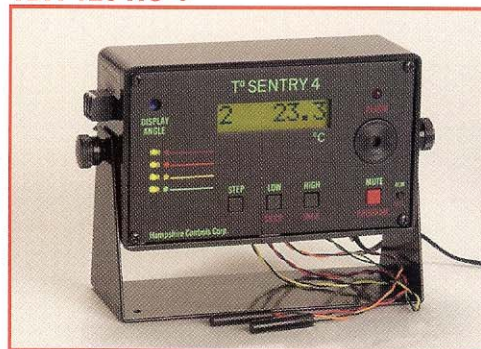

TDA-B

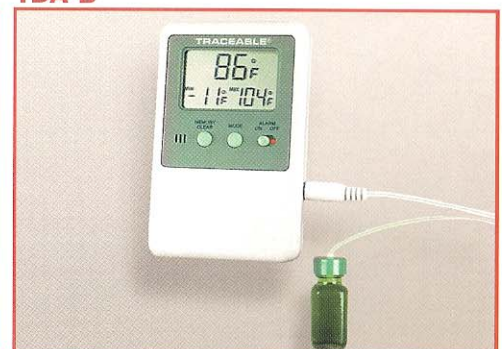

## Temperature Recorders

### FM-4CR120

120 VAC operated flush mounted (factory installed) 4", 7 day circular chart recorder with battery back-up. NIST traceable (not available on all units) 4-CRP: Chart recorder paper 46-CRPEN: Chart recorder pens

### FS-4CRB

Battery operated freestanding self-contained 4", 7 day chart recorder. Unit is portable and can be placed directly in refrigerators and freezers. Perfect for undercounter and chest units. NIST traceable. Powered by a single AA battery. Average battery life of 3 years. Includes pen and one box of 60 charts.


### FS-4CRBD

Same as FS-4CRB but with digital display B4CRP: Chart recorder paper B4-CRPEN: Chart recorder pens


### WM-6CR120

120 VAC operated wall mounted 6", 7 day circular chart recorder. Includes 6' power supply cord and a 15' sensor cable to accommodate temperature probe installation within refrigerator or freezer. Can be surface mounted to nearby walls or work surfaces. Requires simple field installation. 6CRP: Chart recorder paper 46-CRPEN: Chart recorder pens


FM-4CR120

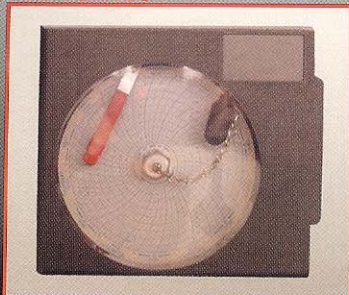

FS-4CRB

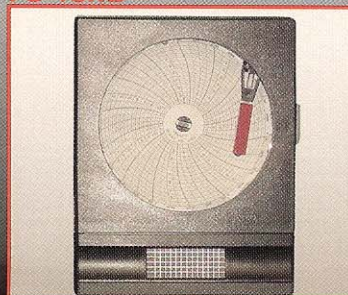

FS-4CRBD

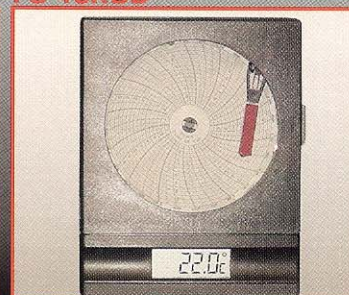

WM-6CR120

