




ENZYME / MANUAL DEFROST FREEZER PLUS SERIES

MANUAL DEFROST & "FAST FREEZE COMPARTMENTS™" ELIMINATE TEMPERATURE SWINGS INSURING ENZYME VIABILITY

Upright Units Featuring "Fast Freeze Compartments™" with inner doors and / or drawers for superior temperature maintenance and recovery. German engineered "Fast Freeze Compartments™" are engineered to utilize inner doors and / or drawers preventing massive cold air loss and resulting product warm-up during exterior door openings. Temperature recovery time is reduced over 80% vs. standard -20° C and -30° C Freezers.

MFP & MCFP Series Standard Features

- Temperature Range: Minimal Avg. Cycle Temperature of -20° C or Lower (MFP-2030 Minimal Average Cycle Temperature of -30° C)
- LR-5 Microprocessor Temperature Controller offers more precise temperature settings and maintenance over mechanical thermostats (see pg. 1)
 - Digital Display
 - Audible and Visual High & Low Temperature Alarm
 - Door Ajar Alarm
 - Battery Back-up for audible and visual alarms
 - Power Failure Alarm
 - Sensor Error Alarm
 - Remote Alarm Contacts
 - High & Low Non-volatile Temperature Memory with reset
 - Audible Alarm Mute Button
- 120 VAC Operation
- Inner "Fast Freeze Compartments™" with doors/drawers for superior temperature maintenance & recovery
- Casters or Caster Platforms
- CFC-Free refrigerant, foam insulation and packaging
- Manual Defrost
- Interior lights (Excludes MFP-0920)
- Door Locks
- Full Length Door Handle
- Magnetic Door Gasket
- Accu-test Certified. Double testing helps ensure accurate temperature performance

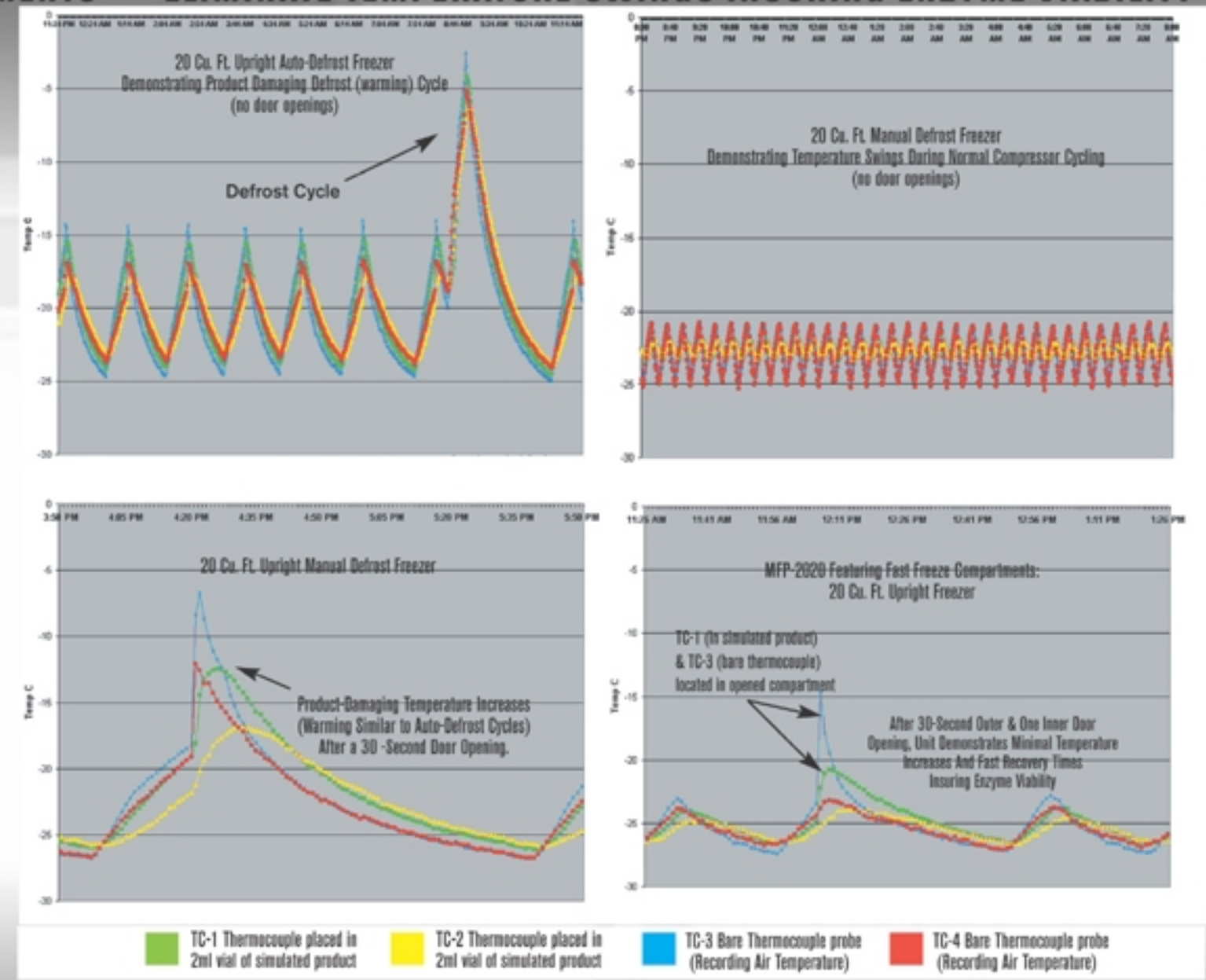

MFP & MCFP Accessories

- (See Accessories Page 20)
- Tower Racks
 - Plastic Storage Boxes & Grids
 - Cardboard Storage Boxes & Grids
 - CryoGloves
 - 120 VAC Operated, Flush Mounted, 4" Circular Chart Temperature Recorder (MFP Units only)
 - Battery Operated Freestanding Internal 4" Circular Chart Temperature Recorder (for chest freezers)



ENZYME / MANUAL DEFROST FREEZER PLUS																						
Model	Cu. Ft.	Doors	Inner Doors/Drawers	Shelves	Usable Int. Dim.			Ext. Dim.			Space between fixed shelves					H.P	Amps	BTU/hr	lbs.	Ship lbs.		
					W"	D"	H"	W"	D"	H"	Top to Bottom											
MFP-0920	9	1 swing	5/2	6 fixed	16	16	50	22	24	6.4*	1.75	9	9.25	6.75	7	7.25	9.5	1/6	7	850	136	176
MFP-2020	20	1 swing	5/0	4 fixed	27	18	57	32	29	7.7*	10.7	10.7	9.5	13.6	9.6			1/4	5	770	315	355
MFP-2030	20	1 swing	5/0	4 fixed	27	18	57	32	29	7.7*	10.7	10.7	9.5	13.6	9.6			1/4	5	770	315	355
MFP-2920	29	1 swing	5/0	4 fixed	32	26	65	35	32	8.4*	11.75	13	12.75	12.8	9			1/3	5	2174	328	368
MCFP-1520	15	1 swing	Chest N/A	2 baskets	42	24	28	48	30	3.5*								1/4	5	770	191	231
MCFP-2020	20	1 swing	Chest N/A	2 baskets	56	24	28	62	30	3.5*								1/4	5	770	220	265